
ogawa
by OGAWA

Revive

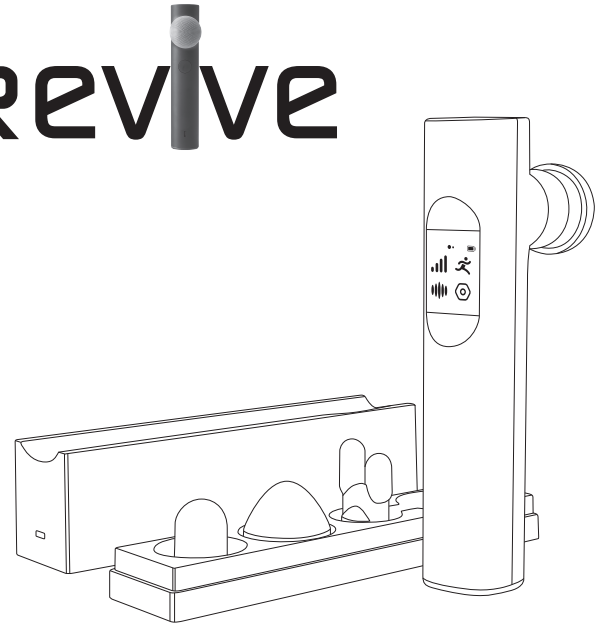User Manual/说明书
OG 1104

*Please provide your order number before contacting customer support

IMPORTANT SAFETY INSTRUCTIONS

WARNING-When using the massager, the following precautions should be followed to avoid damage to the device and to reduce the risk of burns, fire, electric shock or injury to persons:

01. Inspect the massager and read the user manual carefully before use.
02. Do not charge the massager in an environment where temperature is below 0C/32°F or above 40 C/104°F. Charging improperly or at temperatures outside of the specified range may damage the battery and increase the risk of fire.
03. Do not use the massager in an environment where the temperature is higher than 50 C/122°F.
04. Do not expose the massager to fire or excessive temperature. Exposure to fire or temperature above 130°C/266°F may cause explosion.
05. Do not use the massager in an environment where it could be exposed to steam or moisture (such as a bathroom).
06. Under abusive conditions, liquid may be ejected from the battery; avoid contact. If contact accidentally occurs, flush with water. If liquid contacts eyes, seek medical help. Liquid ejected from the battery may cause rrttnburns.
07. Do not use the massager on wet skin.
08. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
09. If the massager is moved too quickly from a low-temperature environment to a warm environment, moisture may accumulate on internal metal parts, which might influence the normal operation of the product, or cause damage to it. To avoid this happening, avoid use for at least one hour after moving it between contrasting temperatures.
10. Do not use a massager that is damaged or modified. A damaged or modified device may exhibit unpredictable behavior resulting in fire, explosion or risk of injury.
11. If you encounter any of the following, please discontinue use of the massager and contact customerservice:
 - a. If it has a damaged charging cord, adapter or charging port.
 - b. If it is not working normally.
 - c. If it has been dropped or damaged.
 - d. If it has been submerged in water.
12. Do not modify or attempt to repair the massager except as indicated in the instructions for use and care.
13. Stop using immediately if you feel discomfort during use.
14. Avoid using on or near the head to avoid injury or health risks.

Menu

Important And Safety Instruction	-----	01
Preparation Before Use	-----	04
Product Diagram	-----	05
Directions For Use	-----	09
Battery Disposal Instructions	-----	14
Cleaning & Maintenance	-----	14
Faulty Inspect And Reset	-----	15
Warranty Claims	-----	16
Contact Us	-----	17

15. Do not use on or near bony prominences.
16. Only use on soft tissue.
17. Stop using if you experience pain or discomfort.
18. Do not use on areas that are hyper-sensitive or have diminished sensitivity due to nerve or circulatory issues.
19. If you fall into any of the following groups, please consult a physician before using the product:
 - a. If you are pregnant or experiencing a menstrual cycle.
 - b. If you have osteoporosis or osteopenia.
 - c. If you have cardiovascular disease.
 - d. If you have neuropathy.
 - e. If you have a pacemaker or other implanted medical devices.
 - f. If you have an acute disease.
 - g. If you have an active/ open skin lesion and/ or trauma to the skin.
 - h. If you are post-operative and still under the direction of a medical practitioner.
20. Avoid using if you are impaired.
21. Close supervision is necessary when this device is used by, on, or near children and physically challenged individuals.
22. Do not fall asleep while operating the massager.
23. Do not use this product while driving.
24. Use this massager only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
25. Do not carry this massager by its charging cord or use its cord as a handle.
26. Do not use outdoors.
27. Prevent unintentional powering of the device. Ensure the switch is in the OFF-position before connecting the device to a power source, picking it up or carrying it. Carrying the device while your finger or fingers are on the ON/OFF switch or connecting the device to a power source with the switch in the ON-position can lead to accidents.

DANGER

To reduce the risk of electric shock:

01. Unplug from the electrical outlet before cleaning.
02. Do not reach for the massager if it falls into water. Unplug immediately.
03. Do not use while bathing or showering.
04. Do not place or store massager where it can fall or be pulled into a bathtub or sink. Do not place in or drop into water or other liquid.
05. Do not use massager in close proximity to loose clothing or jewelry.
06. Keep long hair away from massager while in use.
07. Do not tamper with or alter the unit in any way. There are no user-serviceable components or parts inside.
08. This unit contains a rechargeable Lithium-ion battery, which is neither replaceable nor user-serviceable.
09. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
10. Avoid inserting objects into any opening.
11. To disconnect, power down (press trigger button until it powers down), then remove plug from outlet.
12. Recharge only with the charger specified by the manufacturer. A charger that is suitable for one type of battery pack may create a risk of fire when used with another battery pack.
13. Have servicing performed by a qualified repair person using only identical replacement parts. This will ensure that the safety of the product is maintained.
14. This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance (5V/2A).

PREPARATION BEFORE USE

● Caution before use

Before use, please take off your jewelry to avoid discomfort during massage.
If you feel uncomfortable during initial use or use, please stop using immediately and consult your physician.

● Scope of application


Reduce muscle soreness, promote blood circulation, relieve muscle stiffness, and promote muscle recovery.
This product is only suitable for auxiliary relief, if you need treatment, please go to the hospital for treatment.

● About the status after use

The massage area may turn red for a short time, but it will return to normal after a while.
The massage area may also leave indentation, but after a period of time will return to normal.

When the product is in use, it will stop working or even shut down when it is subjected to excessive resistance. This is the normal motor protection function of the product, and it can be used normally after it is restarted, not because the product is abnormal.
It is normal for the main engine to heat up slightly during charging and use.

PRODUCT DIAGRAM

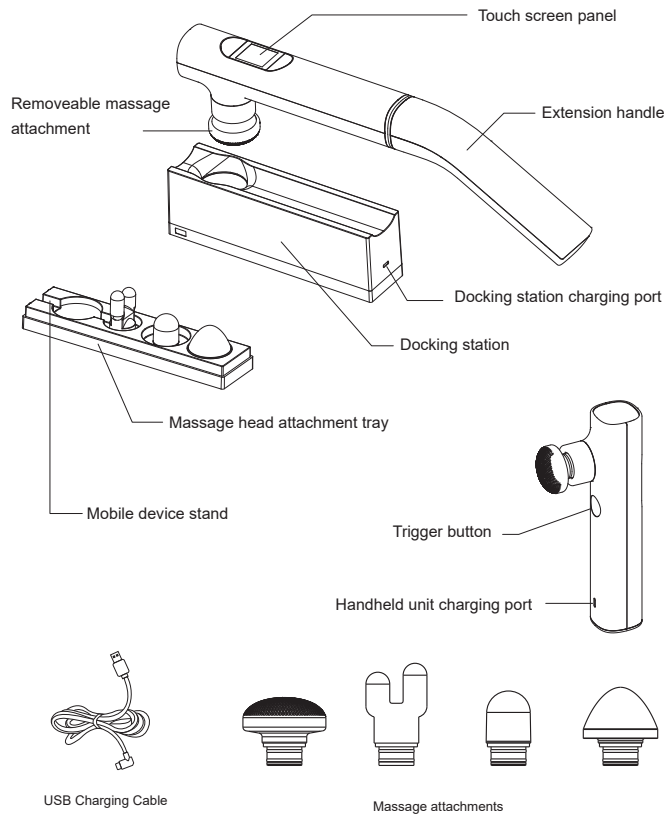PRODUCT DIAGRAM

Dome*1

Suitable for use on most parts of the body.


Cone*1

Soft silicone material, skin friendly soft, suitable for soft tissue muscle use.

Thomb*1

For deep muscle massage such as soles of feet and palms.

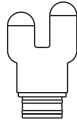Double

Deep muscle massage for large muscle groups.

Note: Bullet massage head and U-shaped massage head material is relatively hard, not suitable for soft tissue massage, in case of muscle damage; When using for the first time, it is recommended to start from the lowest gear, and then raise the gear after the body slowly ADAPTS to the intensity.

Inside the package


Charging Cable

User Manual

Extension Handle

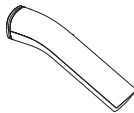
Note: This product does not include adapter, please purchase DC5V 2A adapter.

PRODUCT DIAGRAM


Fascia gun can use wired charging and wireless charging two modes, charging process, charging equipment and products heating is a normal phenomenon.

Wired Charging

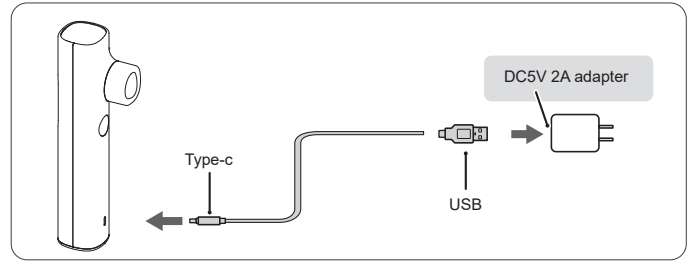Wireless Charging


(Connect the charger to the base and put the product into the base to charge wirelessly)

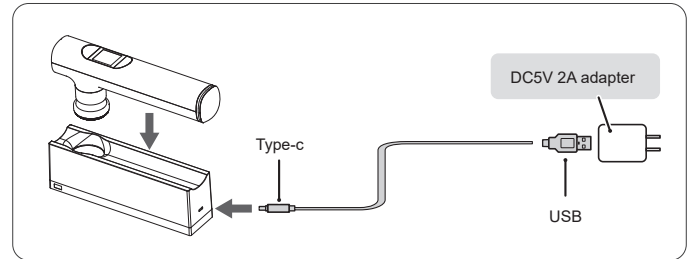The fascia gun cannot be used when charging. During the charging process, the screen will display the charging status. Please remove the charger in time after charging.


PRODUCT DIAGRAM

Before installing the massage head, make sure the product is turned off

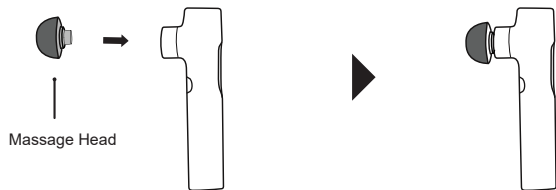
1. Insert the massage head into the mounting hole of the massage head on the fuselage and insert it firmly.
2. When removing the massage head, please hold the root and pull it out straight. If it is more difficult, it can be pulled while rotating, and it can be easier to remove.
3. If the massage head has stains, can use alcohol to wipe clean; Do not wash the massage head or any part of the product.

Mounting extension handle

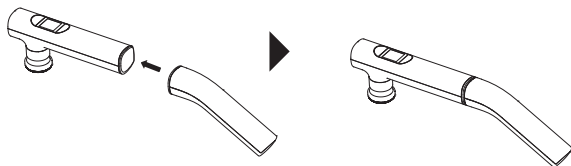Insert the extension handle into the mounting hole at the end of the fuselage and insert it firmly.

DIRECTIONS FOR USE

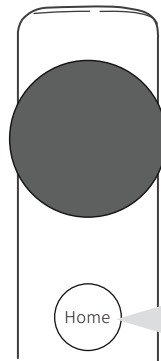Press the Home button to wake up the fascia gun

On state:

1. When the fascia gun is in a non-working state, long press to shut down, short press to return to the upper-level menu;
2. When the fascia gun is in the working state, long press to lock the screen, the fascia gun can continue to work, in the locking screen, short press to wake up the fascia gun;
3. When the fascia gun is in the working state, short press to adjust the massage intensity of the fascia gun.

Note: When the fascia gun is in a non-working state, it can be long pressed to shut down.


The first time the fascia gun is opened, the boot boot will be entered, and the main menu can be entered after the boot setting is completed.

The fascia gun will be rebooted after factory Settings are restored.

⚠ WARNING

- When using for the first time, please press the lowest setting to start using, after the body ADAPTS to the strength, then raise the massage gear; Use only one intensity for scalp and face massage.
- Use only on soft tissues, taking care to avoid nearby bones and joints to avoid damage.
- Due to the inevitable loss of high-frequency impact, it is normal for the product to run at a slightly larger volume. If there is noise or abnormal noise, please contact customer service.

DIRECTIONS FOR USE

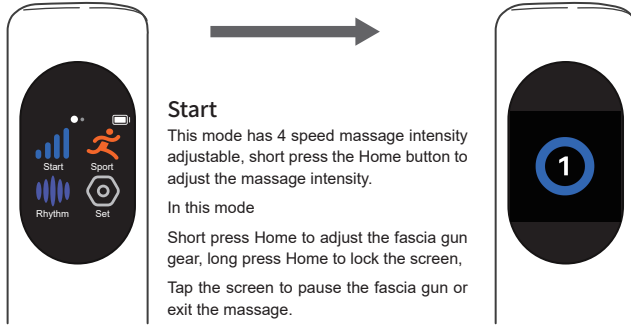Start

This mode has 4 speed message intensity adjustable, short press the Home button to adjust the message intensity.


In this mode

Short press Home to adjust the fascia gun gear, long press Home to lock the screen,

Tap the screen to pause the fascia gun or exit the massage.

Note: In this mode, the product will automatically stop after 10 minutes of operation.

DIRECTIONS FOR USE

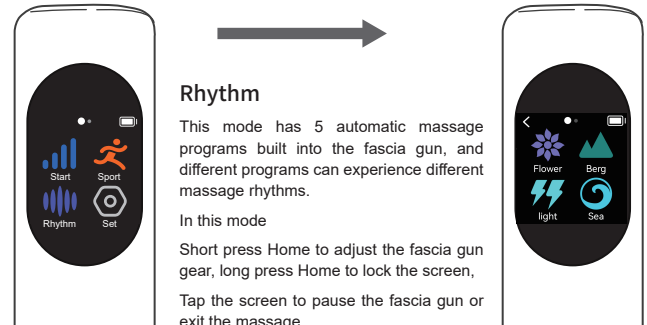Rhythm


This mode has 5 automatic message programs built into the fascia gun, and different programs can experience different message rhythms.

In this mode

Short press Home to adjust the fascia gun gear, long press Home to lock the screen,

Tap the screen to pause the fascia gun or exit the massage.

Note: In this mode, the product will automatically stop after 10 minutes of operation.

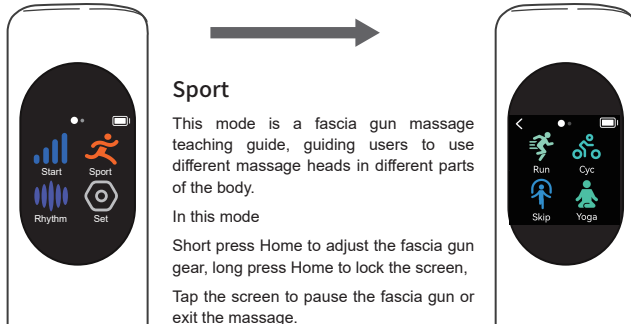
Sport

This mode is a fascia gun message teaching guide, guiding users to use different message heads in different parts of the body.

In this mode

Short press Home to adjust the fascia gun gear, long press Home to lock the screen,

Tap the screen to pause the fascia gun or exit the massage.

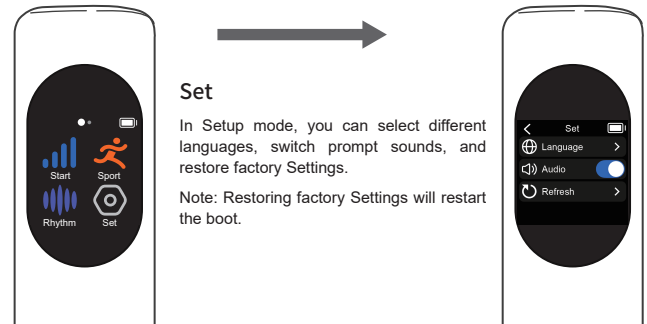
Set

In Setup mode, you can select different languages, switch prompt sounds, and restore factory Settings.

Note: Restoring factory Settings will restart the boot.

DIRECTIONS FOR USE

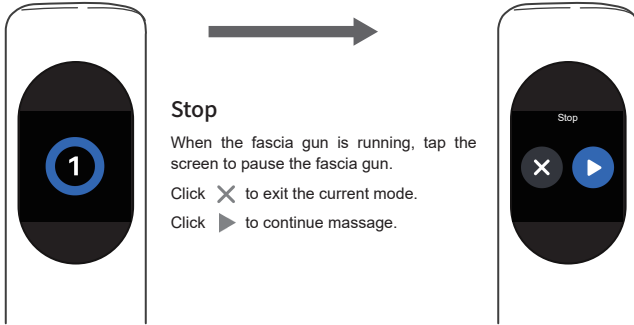
Stop

When the fascia gun is running, tap the screen to pause the fascia gun.

Click  to exit the current mode.

Click  to continue massage.

DIRECTIONS FOR USE

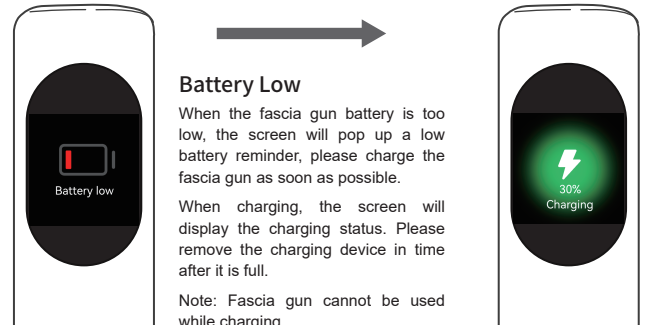
Battery Low

When the fascia gun battery is too low, the screen will pop up a low battery reminder, please charge the fascia gun as soon as possible.

When charging, the screen will display the charging status. Please remove the charging device in time after it is full.

Note: Fascia gun cannot be used while charging.

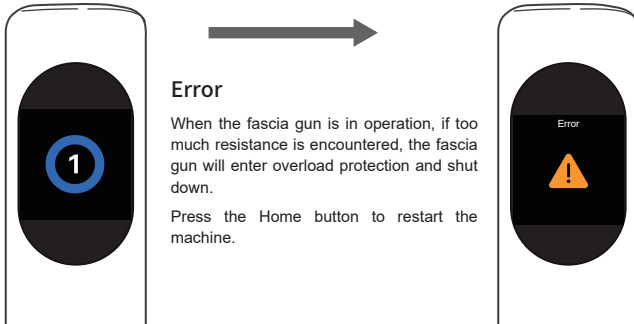Error

When the fascia gun is in operation, if too much resistance is encountered, the fascia gun will enter overload protection and shut down.

Press the Home button to restart the machine.

Product Specifications

Product Name	Revive		
Model No.	OG 1104		
Input Voltage	5V $\overline{\text{---}}$	Power Rating	10W
Current	2A	Battery Capacity	2200mAh
Working Time	10 minutes (product will automatically stop after 10 minutes).		
Product Weight	0.8kg (included all massage head and docking)		
Charging Time	Wired charging, full charge about 4 hour Wireless charging 16H Can use for about 2 hour		

BATTERY DISPOSAL INSTRUCTIONS

The lithium battery in this device is a valuable resource that can be recycled. When the device is at the end of its service life, remove the lithium battery and dispose of it at a facility or establishment that offers battery recycling services. This device should not be discarded with household waste.

CLEANING & MAINTENANCE

Storage

storage over than 3 months, the machine should be recharged after unpacking and inspection. Do not place it close to the floor, walls or roof. Indoor should be well ventilated, avoid intense sunlight illuminate and corrosive gas erosion.

Transport

The product package is equipped with simple anti-shake, suitable for air, railway, road and ship transportation, should avoid rain, snow splashing, inversion and collision.

Power inspection

Check the power supply. If the power supply voltage exceeds the specified range, do not start up and use.

Maintenance

Take measures to prevent moisture, rodents, and pollution. The cleaning should be done in the shutdown state. After use, place the instrument in a safe place, and never put heavy objects on the instrument. When a fault occurs, please ask professional and technical personnel to repair, non-professionals are not allowed to disassemble the product by themselves.

Disposal and Recovery

This product should not be discarded with household waste. When you decide to dispose of this product, please follow the local environmental laws and guidelines.

FAULTY INSPECT AND RESET

Please check if there is a fault in the operation. If the fault cannot be removed or the product still cannot be operated satisfactorily, please cut off the power supply and be repaired by professional technicians. Please do not disassemble the product by yourself!

Basic faults	Possible causes	Basic troubleshooting methods
Not turning on	Battery low	Please charging
Weak massage	Battery low	Please charging

Warranty Claims

1. Warranty is not transferable.
2. This warranty is effective only if the product is purchased from OGAWA.
3. Proof of purchase (original receipt) is required for all warranty repairs.
4. All implied warranties, including but not limited to these implied warranties of fitness and merchantability, are limited to the duration stated in the warranty policy.
5. Damage due to misuse, improper treatment and unauthorized modification and repair are not covered by this warranty.
6. Warranty does not cover accessories and attachments which do not belong to the massage unit.
7. Warranty is not effective to rental, business, commercial, institutional or other nonresidential users.
8. All services covered by this warranty must be approved by OGAWA and repaired by authorized technicians only.
9. If replacement parts for defective materials are not available, OGAWA reserves the right to make substitution in lieu of repair or replacement.

Contact Us

MALAYSIA

Healthy World Lifestyle Sdn Bhd
No.22 Jalan Anggerik Mokara 31/47
Kota Kemuning, 40460 Shah Alam
Selangor Darul Ehsan, Malaysia
Tel: +1800-88-4688
Fax: +603-5121 4386

HONGKONG

OGAWA Health Care
International (HK) Ltd
31/F, YHC Tower,
1 Sheung Yuet Road, Kowloon Bay,
KLN, Hong Kong S.A.R.
Tel: +852 2428 8823
Fax: +85 22410 0119

SINGAPORE

OGAWA Health-Care Pte Ltd
163, Kallang Way #07-11
SINGAPORE 349256
Tel : +065-6846 8588
Fax : +065-6846 7677

VIETNAM

OGAWA Vietnam
171, Nguyen Thi Thap Street
Tan Phu Ward
Dist 7
Ho Chi Minh City, Vietnam
Tel : +84 8 3775 3093 / 3094
Fax : +84 8 3775 3095

PHILIPPINES

O-Healthcare Solution Phil. Inc.
2303, 4th Floor, G&A Building,
Chino Roces Avenue Extension,
Makati City, Philippines
Tel: +632 889 0676/+632 844 6851
Fax: +632 812 0421

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